



Filename: 2026-09-07 - 14.41.22.HSRNA

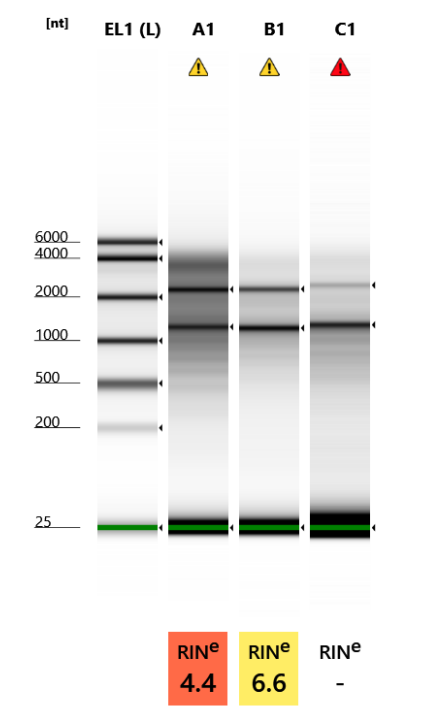

Default image (Contrast 50%), Image is Scaled to Sample

EL1: Electronic Ladder

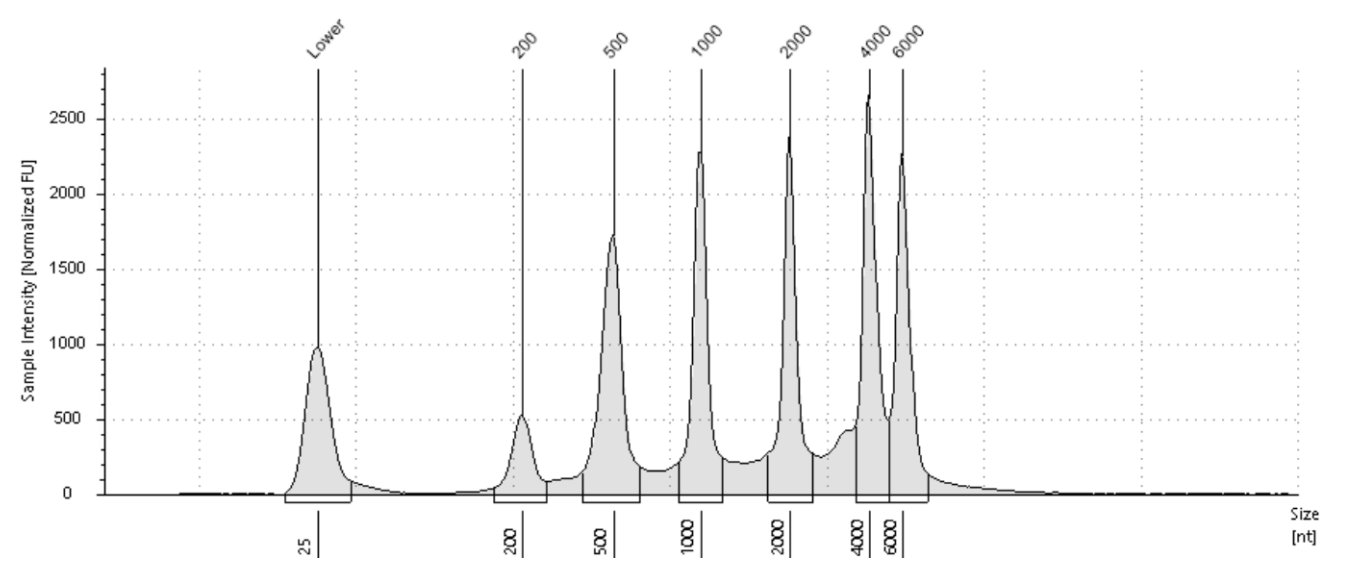


Sample Table

Well	RINe	28S/18S (Area)	Conc. [pg/ul]	Sample Description	Alert	Observations
EL1	-	-	3750	Electronic Ladder		Ladder

Peak Table

Size [nt]	Calibrated Conc. [pg/ul]	Assigned Conc. [pg/ul]	Peak Molarity [pmol/l]	% Integrated Area	Peak Comment	Observations
25	700	700	82400	-		Lower Marker
200	189	-	2770	5.98		
500	638	-	3760	20.25		
1000	557	-	1640	17.66		
2000	580	-	854	18.42		
4000	634	-	466	20.11		
6000	554	-	272	17.58		

A1: POR_R120_H

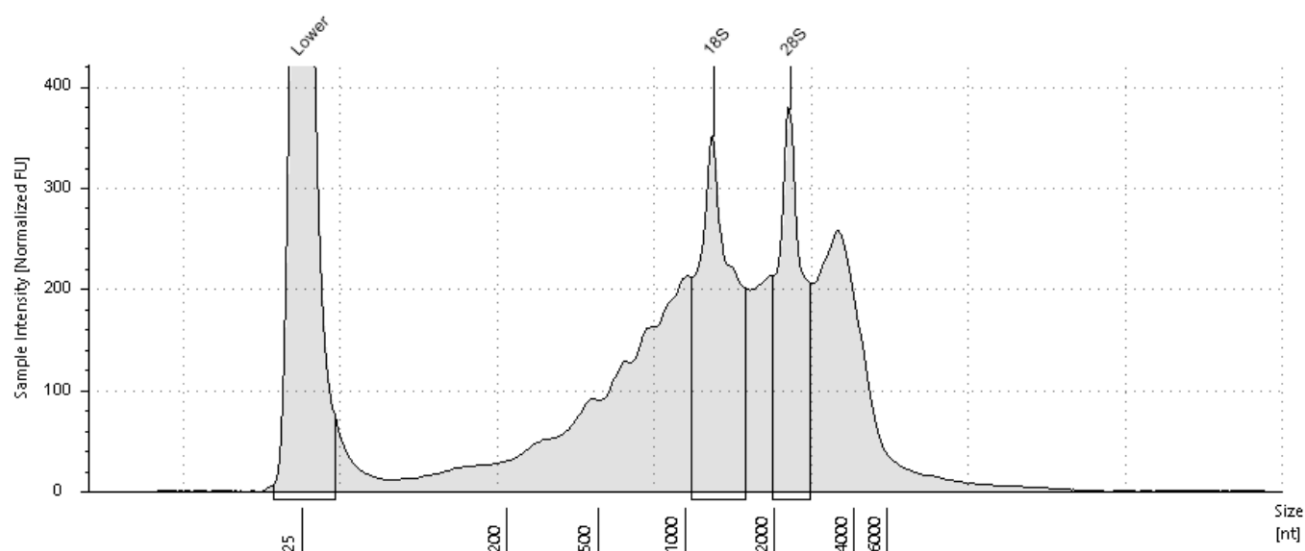


Sample Table

Well	RINe	28S/18S (Area)	Conc. [pg/ul]	Sample Description	Alert	Observations
A1	4.4	0.7	931	POR_R120_H		Caution! Expired ScreenTape device (used after two weeks of first use); RNA concentration outside recommended range for RINe

Peak Table


Size [nt]	Calibrated Conc. [pg/ul]	Assigned Conc. [pg/ul]	Peak Molarity [pmol/l]	% Integrated Area	Peak Comment	Observations
25	700	700	82400	-		Lower Marker
1246	192	-	452	57.61		18S
2299	141	-	181	42.39		28S

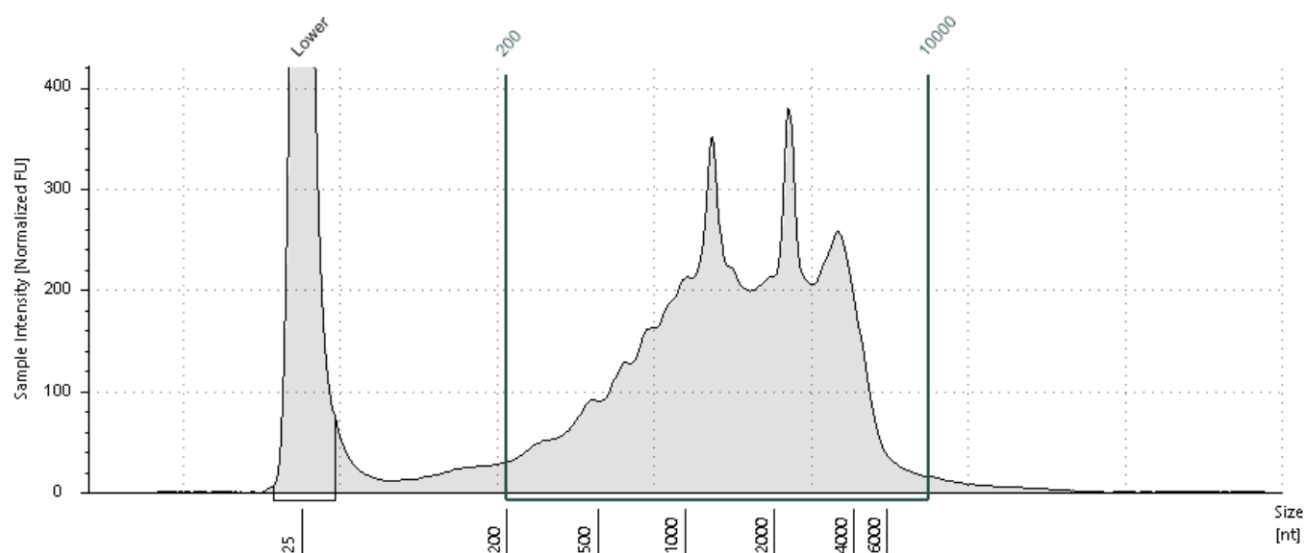

Region Table

From [nt]	To [nt]	Average Size [nt]	Conc. [pg/ul]	Region Molarity [pmol/l]	% of Total	Region Comment	Color
-----------	---------	-------------------	---------------	--------------------------	------------	----------------	-------


200	10000	3188	869	802	93.31	dv200	
-----	-------	------	-----	-----	-------	-------	---

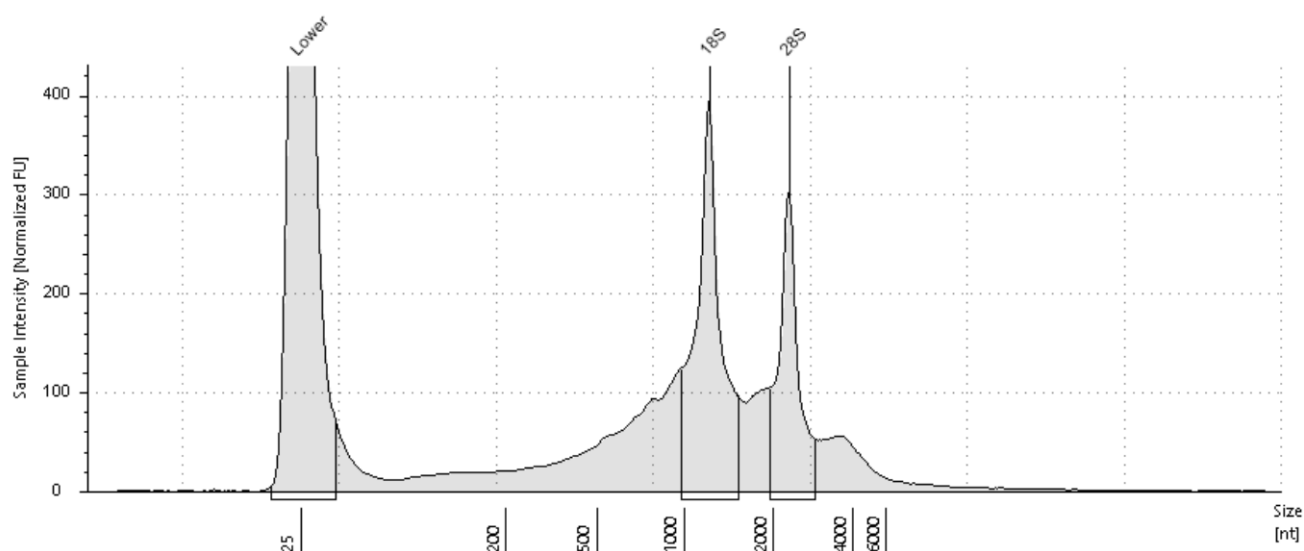B1: POR_R3_C1_Resection**Sample Table**

Well	RINe	28S/18S (Area)	Conc. [pg/ul]	Sample Description	Alert	Observations
B1	6.6	0.6	486	POR_R3_C1_Resection	⚠	Caution! Expired ScreenTape device (used after two weeks of first use); Sample concentration outside functional range for RINe and the assay

Peak Table


Size [nt]	Calibrated Conc. [pg/ul]	Assigned Conc. [pg/ul]	Peak Molarity [pmol/l]	% Integrated Area	Peak Comment	Observations
25	700	700	82400	-		Lower Marker
1220	144	-	348	62.37		18S
2307	87.1	-	111	37.63		28S

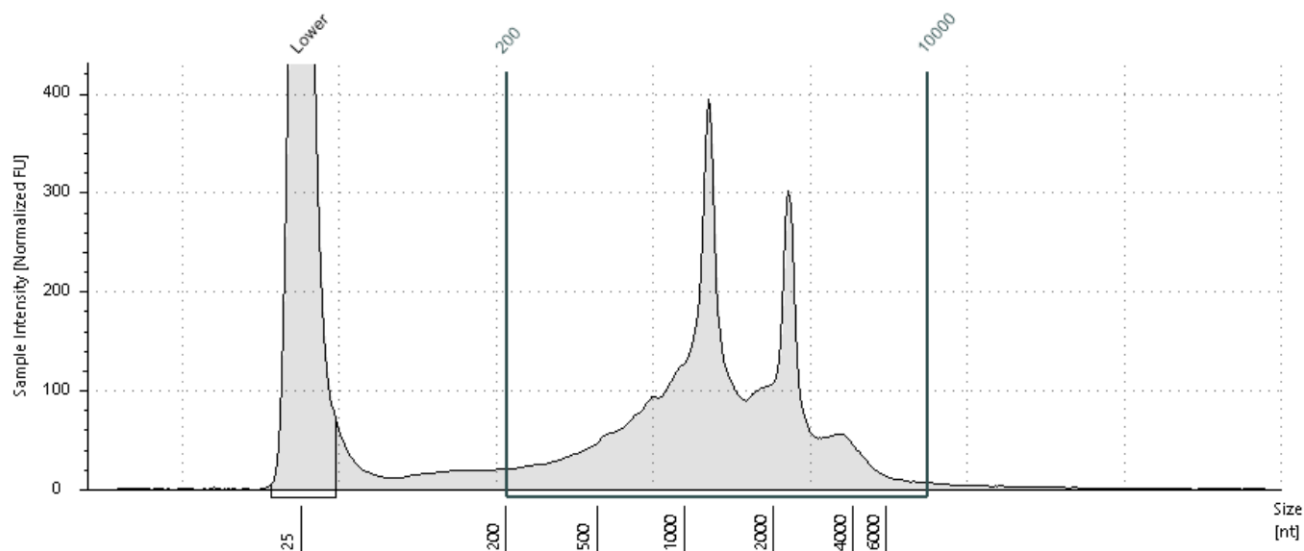**Region Table**

From [nt]	To [nt]	Average Size [nt]	Conc. [pg/ul]	Region Molarity [pmol/l]	% of Total	Region Comment	Color
-----------	---------	-------------------	---------------	--------------------------	------------	----------------	-------

200	10000	2697	436	476	89.73	dv200	
-----	-------	------	-----	-----	-------	-------	---

C1: POR_R3_H3_Resection

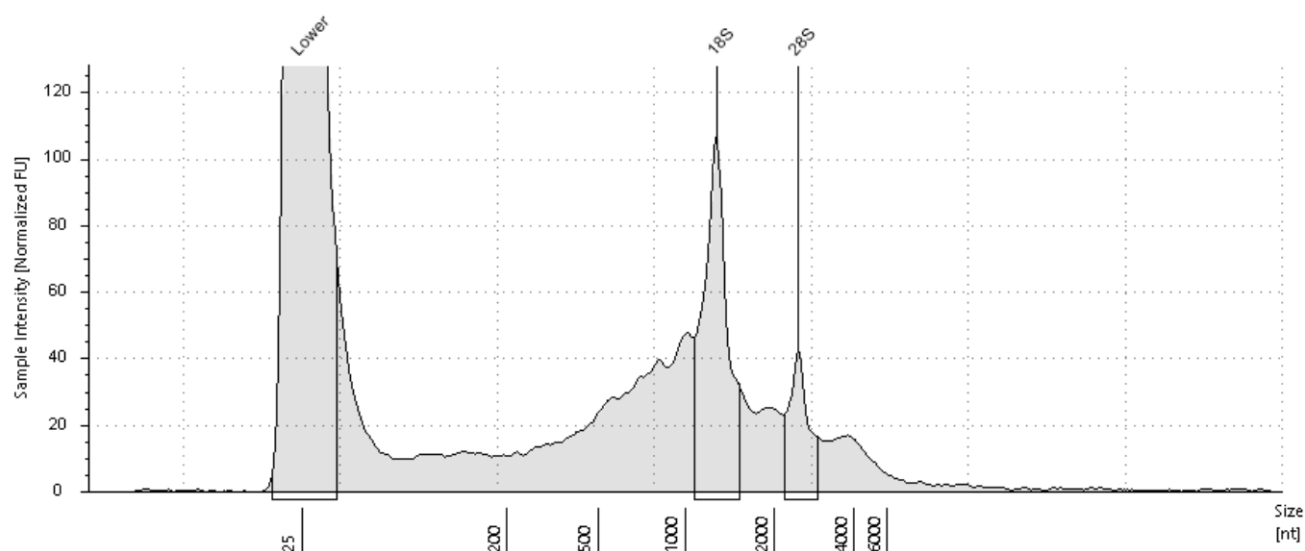


Sample Table

Well	RINe	28S/18S (Area)	Conc. [pg/ul]	Sample Description	Alert	Observations
C1	-	-	171	POR_R3_H3_Resection	▲	Caution! Expired ScreenTape device (used after two weeks of first use); Sample concentration outside functional range for RINe and the assay

Peak Table

Size [nt]	Calibrated Conc. [pg/ul]	Assigned Conc. [pg/ul]	Peak Molarity [pmol/l]	% Integrated Area	Peak Comment	Observations
25	700	700	82400	-		Lower Marker
1286	37.2	-	85.1	76.11		18S
2474	11.7	-	13.9	23.89		28S

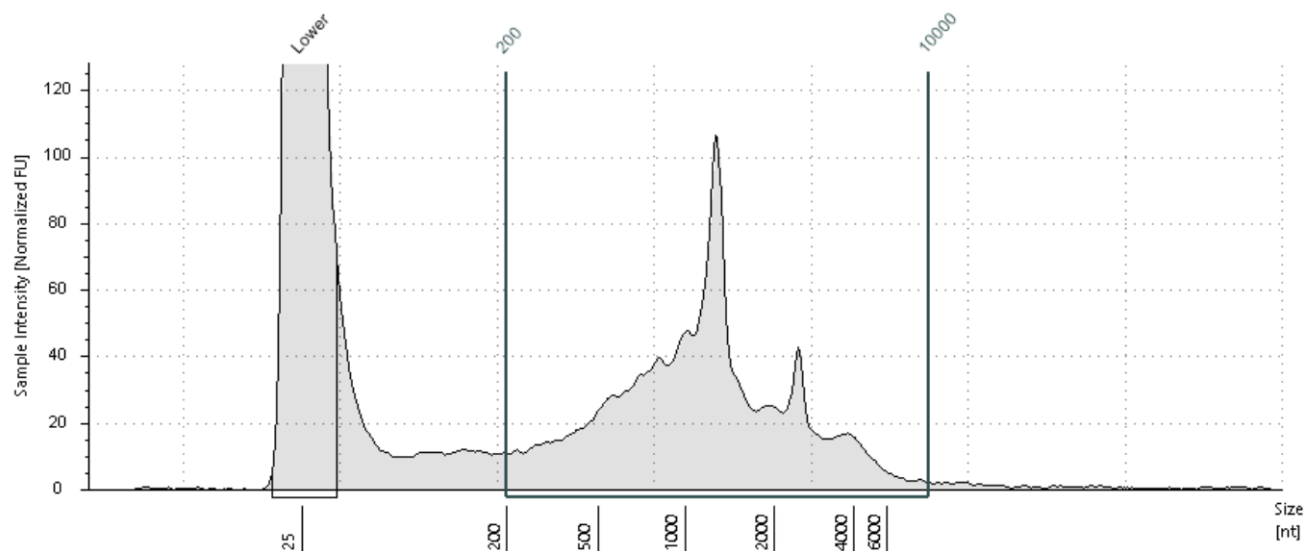

Region Table

From [nt]	To [nt]	Average Size [nt]	Conc. [pg/ul]	Region Molarity [pmol/l]	% of Total	Region Comment	Color
-----------	---------	-------------------	---------------	--------------------------	------------	----------------	-------

200	10000	2787	136	144	79.78	dv200	
-----	-------	------	-----	-----	-------	-------	---

Sample Info

Well	RINe	28S/18S (Area)	Conc. [pg/ul]	Sample Description	Alert	Observations
EL1	-	-	3750	Electronic Ladder		Ladder
A1	4.4	0.7	931	POR_R120_H		Caution! Expired ScreenTape device (used after two weeks of first use); RNA concentration outside recommended range for RINe
B1	6.6	0.6	486	POR_R3_C1_Resection		Caution! Expired ScreenTape device (used after two weeks of first use); Sample concentration outside functional range for RINe and the assay
C1	-	-	171	POR_R3_H3_Resection		Caution! Expired ScreenTape device (used after two weeks of first use); Sample concentration outside functional range for RINe and the assay